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peqlab

CONTENT

WARRANTY	1
SAFETY PRECAUTIONS	1
PACKAGING LIST	2
INSTALLATION	2
OPERATING ELEMENTS	3
OPERATING FEATURES	3
CALIBRATION	4
CLEANING AND MAINTENANCE	5
TECHNICAL DATA	5
TROUBLESHOOTING	6
TECHNICAL SUPPORT AND ORDERING INFORMATION	7
DISPOSAL INSTRUCTION	7

WARRANTY

PEQLAB guarantees that the system you have received has been thoroughly tested and meets its published specification.

However, immediately upon arrival, please check carefully that the shipment is complete and has not been damaged in transit. For missing parts or to report any kind of damage, please contact PEQLAB (see 'TECHNICAL SUPPORT AND ORDERING INFORMATION'). Please retain all packaging materials until the delivery has been completely checked since this will speed up the return of goods if required and reduce environmental impact. Any form of returns, replacements or credit notes must be agreed in advance by PEQLAB.

For the PerfectBlot Hybridisation Oven PEQLAB guarantees a warranty period of 24 months if the products have been used solely according to the instruction manual and if not agreed differently. No liability is accepted for loss or damage arising from incorrect use. PEQLAB's liability is limited to the repair or replacement of the unit or refund of the purchase price, at PEQLAB's discretion. PEQLAB is not liable for any consequential damages.

PEQLAB reserves the right to alter the technical specifications without prior notice. This will enable us to implement developments as soon as they arise.

SAFETY PRECAUTIONS

1. Although the PerfectBlot Hybridisation system has been designed to minimize user exposure to radioactive materials, proper precautions must be taken when using radioactive material.
2. **Always inspect the hybridisation bottles, caps and seals before use.** Do not use bottles, caps or seals that are cracked or chipped. Cap seals should be replaced if they are severely deformed or discolored. Special attention should be paid to inspecting bottle rims and threads.
3. **DO NOT take bottles to temperatures higher than 70 °C without opening and retightening the cap to relieve pressure build-up.** Temperatures above 70 °C can cause breakage of hybridisation bottles. If chamber temperatures rise above 70 °C when bottles are in use, **DO NOT OPEN THE OVEN DOOR.** Reduce the temperature setting to below 70 °C or turn the power off and allow sufficient time for the bottles to cool below 70 °C before opening the door.

PACKAGING LIST

The package should include:

- PerfectBlot Hybridisation Oven
- Rotisserie
- 4 adjustable feet
- Drip tray
- Rotisserie brass locking pin
- 2 hybridisation bottles (300 x 35 mm)
- 1 package of mesh
- Manual

INSTALLATION

Remove the drip tray, adjustable feet and rotisserie from the chamber (refer to OPERATING FEATURES for Rotisserie Removal and Installation).

Install the four adjustable feet and locate the Hybridisation Oven on a stable, flat surface near a grounded electrical outlet. The location selected should be out of direct sunlight and away from heat producing sources or hot or cold air drafts. At least 5 cm ventilation clearance is required around all sides of the oven. Level the incubator using the four adjustable feet. Clockwise rotation of a foot raises the oven. You may wish to use a level placed in the chamber for optimum leveling. Check the rotisserie and set it for horizontal or angled use by adjusting the screws in the rotisserie hub (see OPERATING FEATURES, Setting Rotisserie Offset Angle). Plug the unit in to a properly rated and grounded electrical outlet. It is now ready for use.

To stack ovens: The PerfectBlot Ovens can be stacked two high. When stacking, it is advisable to use either rubber furniture cups under the feet of the upper unit or flat rubber matting at each corner of the upper unit or across the entire upper surface. You may wish to remove the adjustable feet from the upper unit when stacking. The adjustable feet can slide off the top of the lower unit and/or scratching can occur if rubber cups or matting are not used.

To stack a unit, place the upper oven on top of the lower oven slightly toward the back such that the upper oven does not prevent the lower oven door from moving freely.

OPERATING ELEMENTS

The controls for PerfectBlot Ovens include power to the oven, temperature setting and display and rotisserie power and speed. The controls are located at the bottom front of the oven. A safety electrical circuit breaker is located on the back of the oven.

Main Power Switch

This switch (illuminated green when on) turns the power to the unit off and on.

Temperature Controller and Temperature Set

The controller has a 3-digit display for displaying chamber temperature or set point information. UP and DOWN arrow pads are used to change the set point and controller mode of operation. To enter the set point mode of operation, press either the UP or DOWN arrow pad once. The display will start to blink, going from bright to dim. While blinking, the display is showing the set point. To change the set point, use the UP and DOWN arrow pads. If the arrow pads are not pressed for 5 seconds, the display will stop blinking and will read the chamber temperature. Allow at least one hour for the chamber temperature to stabilize and 24 hours for optimum stabilization.

Heating Indicator

This indicator will illuminate green when the controller is calling for heat from the heater. This indicator will be on continuously while the oven heats up to the set temperature and will then cycle on and off at the set temperature.

Rotisserie Power Switch

This switch turns on the power to the rotisserie and is used to start and stop the rotisserie.

Rotisserie Speed

This knob controls the rotisserie speed. Clockwise rotation increases speed. Nominal speed range is 4 to 20 rpm.

OPERATING FEATURES

Rotisserie

The stainless steel rotisserie supplied with the unit will hold either twelve large (300 mm) bottles or twenty four small (150 mm) bottles. The clips that hold the bottles may be squeezed inward to accommodate 50 ml tubes. Optional rotisseries that hold other styles of bottles or vertical 50 ml tubes are also available. Always load the rotisserie with an even number of bottles in a balanced pattern.

Setting Rotisserie Offset Angle

The rotisserie can be adjusted to hold bottles either horizontally or at a slight angle for more vigorous wetting action. To adjust the angle, first loosen the screws on one hub then rotate the hub to create an angle offset between the two hubs, then retighten the screws. Additional angle can be obtained by loosening the screws and rotating the second hub in the opposite direction from the first hub.

Rotisserie Removal and Installation

Caution! Rotisserie and oven sides may be hot! The rotisserie is held and driven by a two-prong fork drive on one side and locked in a cradle on the other side by a brass locking pin. To remove the rotisserie, first remove the brass locking pin, then lift the rotisserie up from the cradle while simultaneously moving the rotisserie away from and off the two-prong fork drive. Rotisserie installation is the reverse. Bottles may need to be removed before attempting to remove or install the rotisserie.

Drip Tray

Each oven comes with a drip tray that is removable for easy clean up of spills.

Rocking Platform Installation

An optional rocking platform is available for use in the PerfectBlot oven chamber. Before installing the rocking platform, the rotisserie must first be removed (see above). To install the rocking platform, first install the fork adapter (rectangular bar with holes) on the rotisserie drive fork by pressing the adapter over the two prongs and tighten the Allen screw in the end of the adapter using the wrench provided. Do not over tighten.

Next, use the rotisserie speed control to bring the fork drive and adapter to a horizontal position. Place the rocking platform in the oven and connect the platform drive arm to the fork drive adapter using the thumbscrew. Two positions are available on the adapter. Select the position that allows the platform to sit horizontal (or close to horizontal) when the fork drive and adapter are horizontal.

The platform drive arm can be adjusted and the second hole in the adapter used to obtain more or less rocking angle or to achieve a perfectly horizontal platform at mid-rock. If the length of the arm is adjusted, be sure to retighten the locking nuts on the arm.

To operate the rocking platform, turn the rotisserie power switch to the ON position and adjust the speed using the rotisserie speed control knob.

CALIBRATION

The oven is factory calibrated at 65 °C. The oven can be re-calibrated after the chamber temperature has stabilized at the set point for several hours. Suspend a certified reference thermometer in a hybridisation bottle such that it does not touch the glass bottle sides or (less accurate) tape the reference thermometer to the outside wall of a hybridisation bottle. Mount the bottle on the rotisserie (and a second bottle for balance) and run the unit at the desired temperature for 2 hours. Compare the oven display to the reference thermometer. If there is an unacceptable difference, place the controller into calibration mode by pressing both the UP and DOWN arrow pads at the same time until the two outside decimal points begin to flash. While the decimal points are flashing, the display can be calibrated to match the reference thermometer by pressing the UP or DOWN arrow pads until the display reads the correct value. Allow the incubator to stabilize again, and recalibrate if necessary.

CLEANING AND MAINTENANCE

No routine maintenance is required for the electrical or mechanical components of the unit. The incubator exterior, interior, rotisserie and rocking platform (if present) should be wiped down periodically with a soft damp cloth with mild soap. Do not use chlorine-based bleach or abrasives. Any spills in the incubator and/or on the rotisserie, drip pan or rocking platform should be cleaned immediately. **Be sure to disconnect the power cord before cleaning or decontaminating the oven.**

Decontamination

Should a spill occur in the oven, it is easily cleaned by wiping the affected area first with a dilute solution of IsoClean or CountOff (or similar decontamination agent) and then with distilled water. This method works on the oven interior, drip tray, rocking platform and rotisserie. Bottles and caps may also be soaked in a dilute decontaminating detergent such as IsoClean or CountOff and then rinsed in distilled water.

TECHNICAL DATA

Temperature range	5 °C over ambient to 80 °C
Temperature accuracy / uniformity	± 0.1 °C / ± 0.5 °C
Capacity	12 long or 24 short hybridisation bottles (300 x 35 mm, 150 x 35 mm)
Exterior (L x H x D)	483 x 495 x 445 mm
Chamber (L x H x D)	368 x 343 x 279 mm
Weight	26.8 kg
Electrical	230 V, 50 - 60 Hz, 3 amps 120 V, 50-60 Hz, 6 amps

TROUBLESHOOTING

Below are some possible solutions to potential problems. If these suggestions are unsuccessful, do not hesitate to contact us (see 'TECHNICAL SUPPORT AND ORDERING INFORMATIONS').

Problem	Solution
Oven will not power up or not heat	Check power cord, outlet and unit circuit breaker.
Rotisserie will not turn	Be sure the rotisserie power is turned on, increase the speed.
Rotisserie speed erratic, jumps	Check for balanced load - an even number of bottles should be loaded in opposing positions.
Temperature too high	Check set point and calibration. Readjust if necessary.
Chamber temperature rises above set	During initial heating in normal operation, the chamber temperature will rise above set point and then settle back down. If the door is opened for a long period of time, this may also happen.
Temperature will not remain stable or display shows "LO"	Check that set point is at least 5 °C above ambient.
Indicated temperature is unstable	A slight variation of +/-0.1 °C is normal. Larger fluctuations may be ambient variations from drafts, door opening and closing, a fan obstruction or failure or electrical interference from RFI (motors, compressors, etc.).
Temperature is too low	If door has been opened, unit may not have recovered yet. Confirm temperature set point.
Unit will not heat above a certain temperature	Confirm set point. Confirm with a thermometer that actual chamber temperature is below set point. Recalibrate if necessary.
Temperature display and reference thermometer do not match	Be sure that unit has been allowed to stabilize for 1 hour. Thermometer should be suspended in a rotating bottle and not touching bottle surface. Only certified reference thermometers should be used.
Can not adjust set point or calibration	Turn unit off for 5 seconds to reset. If problem persists call service.
Unit not calibrated at the temperature and load most often used	For best accuracy, temperature should be recalibrated at the temperature and load most often used.

TECHNICAL SUPPORT AND ORDERING INFORMATION


For technical questions and more detailed information on PEQLAB's products please visit www.peqlab.com to find the respective contact person.

Please have the unit's serial number (located on the back panel of the instrument) available when calling. Should an item require return to PEQLAB for service, a decontamination form must be completed first by the user. Items sent without decontamination form will not be accepted.

All returns must be pre-approved by PEQLAB!

Item	Description	Cat. No.
PerfectBlot Hybridisation Oven, 230 V	Capacity for 12 long or 24 short hybridisation bottles	91-5510
PerfectBlot Hybridisation Oven, 110 V	Capacity for 12 long or 24 short hybridisation bottles	91-5510-US
Rotisserie	for 16 x 50 ml tubes in vertical position	91-5511
	for 4 hybridisation bottles (300 x 70 mm)	91-5513
Rocking platform	Rocking frequency 8 - 40/min	91-5512
Hybridisation bottle	150 x 35 mm incl. cap	91-5112
	300 x 35 mm incl. cap	91-5113
	300 x 70 mm incl. cap	91-5114

DISPOSAL INSTRUCTION

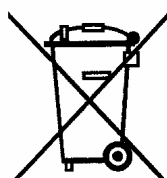This product shouldn't be disposed of with regular waste. Please contact PEQLAB or your local distributor for disposal of this instrument. Waste electronic equipment is not useless waste. Environmentally friendly disposal can help to retrieve valuable resources. Thereby you contribute to the protection of the environment and human health.

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D PEQLAB Biotechnologie GmbH, 91052 Erlangen, Freecall (D): 0800 100 20 16, info@peqlab.de, www.peqlab.de
AT PEQLAB Biotechnologie GmbH, 6404 Polling, Tel: +43 (0) 5238 84 169, info@peqlab.at, www.peqlab.at
UK PEQLAB Ltd., Southampton SO31 7ZN, Freecall (UK): 0808 202 1302, info@peqlab.co.uk, www.peqlab.co.uk
USA PEQLAB LLC, Wilmington, DE 19810, Toll-Free (US): 877 737 5220, info@peqlab.us, www.peqlab.us